
ΔVΔNOS*




CORFLO* NASOGASTRIC / NASOINTESTINAL FEEDING TUBES

NG/NI
Tubes



Instructions for Use



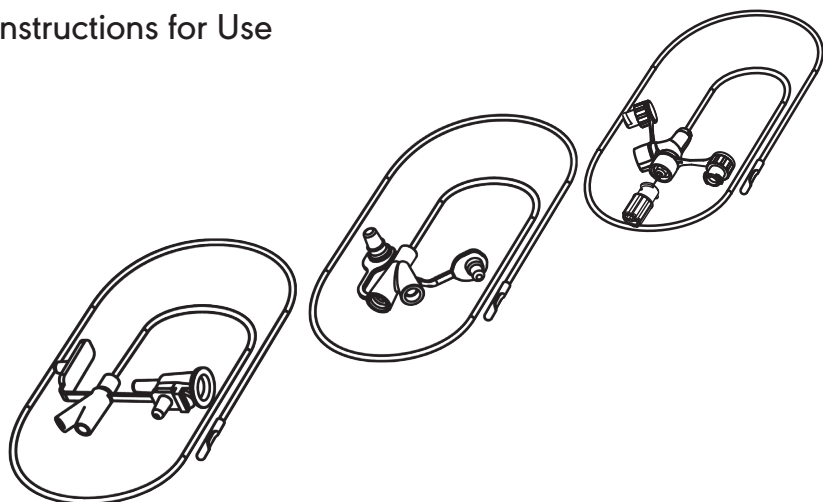
Universal Connector

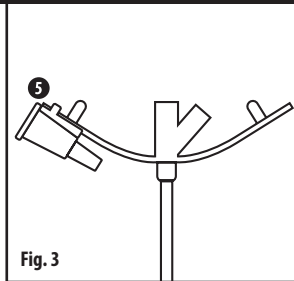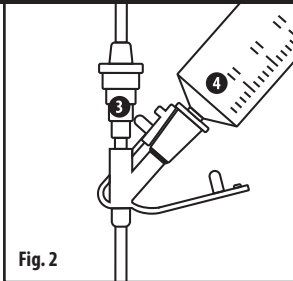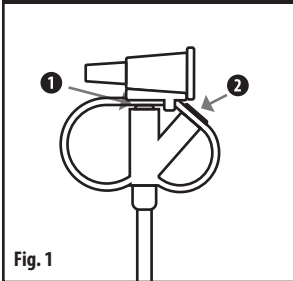ANTI-IV® Connector

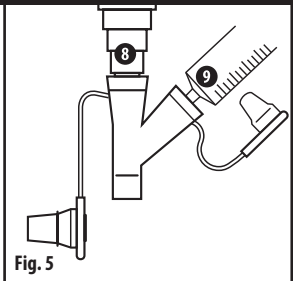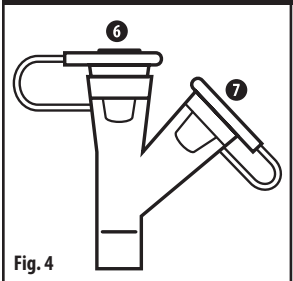ENFit® Connector

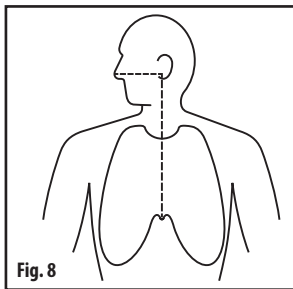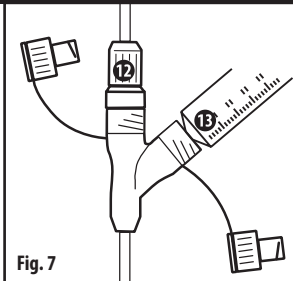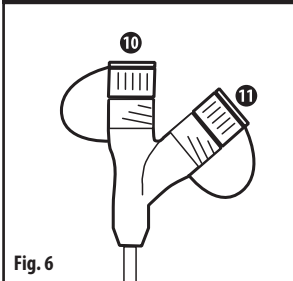Diameter	Length	Single Use Only
Do not use if package damaged	Not made with natural rubber latex	Product is not made with DEHP as a plasticizer
Product is not made with BPA	Product is not made with PVC	Rx Only
MR Conditional	Caution	Consult instructions for use

AVANOS* CORFLO*

Nasogastric / Nasointestinal Feeding Tubes

Rx Only: Federal Law (USA) restricts this device to sale by or on the order of a physician.

For use by Healthcare Professionals only.

Description

The AVANOS* CORFLO* Nasogastric/Nasointestinal Feeding Tube is a medical grade polyurethane feeding tube that has been specifically designed for patient comfort and safety during tube insertion and use. It is intended for use in patients who require intermittent or continuous tube feedings via the nasogastric or nasointestinal pathway. Options include sterile/non-sterile, weighted/non-weighted, and with/without stylet.

Available Connectors

- Universal Connector (Accepts Luer and Oral syringes)
- ANTI-IV* Connector
- ENFit* Connector

Indications for Use

The AVANOS* CORFLO* Feeding Tube is intended for use in those patients who require intermittent or continuous tube feedings via the nasogastric or nasointestinal route.

Device Lifetime

Device has been validated for use up to 30 days.

Warnings:

- ⚠ **The patient should not lean forward, nor should the head and neck be extended.**
- ⚠ **Premeasurement of tubing length is essential. Do not insert excess. Occlusion may result from kinking of tube.**
- ⚠ **Coughing may indicate passage of tube into trachea. If tracheal passage is suspected, remove tube. Absence of coughing does not confirm tube placement in the stomach. If resistance is encountered, immediately remove tube. Notify clinician. Care should be taken if any type of endotracheal device is in place, as it may guide feeding tube into trachea. Misplacement of the feeding tube into trachea or lungs may result in serious injury.**
- ⚠ **Tube position in the stomach must be confirmed prior to flushing and use.**
- ⚠ **Never reinsert stylet when tube is in patient.**
- ⚠ **Vigorous syringe force should not be used to irrigate, administer liquids or unblock the tube.**

Tube Insertion

This feeding tube is to be inserted by trained and competent individuals or clinicians, following institution/facility/hospital protocols.

1. Explain procedure to patient if applicable.
2. Position patient in sitting or Fowler's position as tolerated.
 - ⚠ **Warning: The patient should not lean forward, nor should the head and neck be extended.**
 - ⚠ **Caution:** Do not use a stylet if the tube is not packaged with one, as the stylet may exit the distal tip.
3. Remove tube from package. Close access ports.
4. Measure length of tube to be inserted to ensure that tip/bolus enters the gastric region. Place exit port of tube at tip of nose. Extend tube to earlobe, then to xiphoid process (**Fig 8**). Use the printed centimeter marks on the tube to aid intubation and check for tube migration.
 - ⚠ **Warning: Premeasurement of tubing length is essential. Do not insert excess. Occlusion may result from kinking of tube.**
5. Determine preferred nostril for insertion. Provide cooperative patient with glass of water and straw.

6. Prior to insertion, activate lubricant on tip by dipping in water.
7. Direct tube posteriorly, aiming tip parallel to nasal septum and superior surface of hard palate. Advance tube to nasopharynx, allowing tip to seek its own passage. As patient swallows sips of water, advance tube through esophagus into stomach with gentle motion.
8. Utilize printed centimeter marks on the tubes to aid insertion.
 - ⚠ **Warning: Coughing may indicate passage of tube into trachea. If tracheal passage is suspected, remove tube. Absence of coughing does not confirm tube placement in the stomach. If resistance is encountered, immediately remove tube. Notify clinician. Care should be taken if any type of endotracheal device is in place, as it may guide feeding tube into trachea. Misplacement of the feeding tube into trachea or lungs may result in serious injury.**
9. Confirm tube position and integrity per institution/facility/hospital protocol (e.g. pH measurement, X-ray, etc.).
 - ⚠ **Warning: Tube position in the stomach must be confirmed prior to flushing and use.**
10. For tubes packaged with a stylet: After tube position in the stomach is confirmed, remove stylet by flushing tube through the side port with up to 10 ml of water to activate internal lubricant immediately prior to stylet removal.
 - ⚠ **Warning: Never reinsert stylet when tube is in patient.**
11. Secure tube per institution protocol.

Tube Maintenance

1. Follow your institution/facility/hospital protocol or clinician's order.
2. It is recommended the tube be irrigated every 4 hours with up to 20 ml of water (up to 10 ml for infants or children) before and after medication administration or when feeding formula is interrupted.
 - ⚠ **Warning: Vigorous syringe force should not be used to irrigate, administer liquids or unblock the tube.**
3. The feeding tube should be monitored, regularly assessed, and replaced when clinically indicated based on functionality and patient condition.

To Unclog a Tube

1. Do not use excessive force to flush the tube. Excessive force can perforate the tube and can cause injury to the gastrointestinal tract or risk of aspiration.
2. If the tube becomes clogged, it is best to first attempt clearance immediately by using aspiration with a sterile syringe. Use a 30 to 60 ml syringe. Do not use smaller size syringes as this can increase pressure inside the tube and potentially rupture the tube.
3. If aspiration fails, you may attempt gentle flushing. Use a 30 to 60 ml syringe with sterile fluid. Do not use air, as it is compressible.
4. If gentle flushing is not sufficient to clear the clog, consider use of an enzymatic-based declogging solution (e.g. CLOGZAPPER®).
5. If the tube cannot be unclogged easily (without force), replace the tube with a new one.

Note: When administering crushed medicine, be sure to crush it finely and dilute it as appropriate (until there are no visible particles) to ensure the solution is of a low enough viscosity to avoid clogging.

AVANOS* CORFLO* Nasogastric / Nasointestinal Feeding Tubes with Universal Connectors

- **Tube Capped (Fig. 1)**
 - 1 Feed Port
 - 2 Auxiliary Port

- **Tube Connected (Fig. 2)**
 - ③ Administration Set
 - ④ Syringe
- **Tube Uncapped (Fig. 3)**
 - ⑤ Irrigation Syringe Adaptor

Features

- Medical grade polyurethane
- Product identification printed on tube
- Centimeter markings (approx.) printed on tube to aid placement and check migration
- Radiopaque tube and tip
- Tungsten filled insertion guide tip (weighted tubes only)
- Water-activated C-19* lubricant on tip
- Available with an exclusive anti-clog exit port
- Separate access port for irrigation, aspiration, or additives
- Irrigation adapter color coded for length

AVANOS* CORFLO* Nasogastric / Nasointestinal Feeding Tubes with ANTI-IV* Connectors

- **Tube Capped (Fig. 4)**
 - ⑥ Feed Port
 - ⑦ Auxiliary Port
- **Tube Connected (Fig.5)**
 - ⑧ Administration Set
 - ⑨ Syringe

Features

- ANTI-IV* Y connector on proximal end of feeding tube will not accept Luer Lock or standard Luer slip syringes
- Product identification printed on tube
- Centimeter markings (approx.) printed on tube to aid placement and check migration
- Radiopaque tube and tip
- Tungsten filled insertion guide tip (weighted tubes only)
- Water-activated lubricant on tip
- Exit ports:
 - Anti-clog bolus
 - Bolus
- Side port for irrigation, aspiration, or additives
- Color coded for length

Note: When withdrawing aspirate, use a catheter tip syringe in main port or an oral syringe in side port and slowly withdraw plunger.

AVANOS* CORFLO* Nasogastric / Nasointestinal Feeding Tubes with ENFit® Connectors

- **Tube Capped (Fig. 6)**
 - ⑩ Feed Port
 - ⑪ Auxiliary Port
- **Tube Connected (Fig. 7)**
 - ⑫ Administration Set
 - ⑬ Syringe

Features:

- Distinctive YELLOW* tube design
- Exit ports:
 - Anti-clog bolus
 - Bolus
- ENFit® design employed for reducing the risk of tubing misconnections. ENFit® satisfies the ISO Standard (80369-3) for small bore enteral devices.
- Water-activated C-19* lubricant on tube tip
- Dual port to maintain closed system
- For YELLOW* feeding tubes Radiopacity: 20% barium tube, 40% barium tip

Ability to Connect to Non-Enteral Medical Devices

- The ENFit® connectors were designed to prevent misconnections between enteral devices and other devices used in various medical applications. However, the design of the ENFit® connector cannot overcome all chances of misconnection. The following connector types are potential misconnections for the ENFit® connector (feeding/ medication access port) of this enteral feeding tube:
 - Suction ports on Endotracheal Suction Systems
 - Respiratory circuit filtration connectors
 - Oxygen inlet connectors for Resuscitation Devices
 - Baxter IV Solution Bag ports (such as NaCl, Ringers Solution, etc.)
 - Sample ports on drainage bags
 - Peritoneal Dialysis connectors
 - Cones & sockets of ISO 5356-1:2004 & ISO 5356-2:2004
 - Temperature sensor connectors & mating ports of ISO 8185:2007
 - Oxygen nipples as defined in EN 13544-2:2002

CORFLO® NG/NI feeding tubes have satisfied biocompatibility testing as a device for long-term use per ISO 10993-1.

Note: Placement and use of any feeding tube may result in patient discomfort.

Complications

Infrequently occurring factors as a result of use or misuse of any feeding tube may include:

- Pneumothorax
- Aspiration
- Tissue irritation or necrosis
- Contamination
- Delay or misuse of medication or nutrition and related complications or the need for additional medical procedures
- GI perforation
- Airway obstruction
- Allergic reaction
- Delay in diagnosis

MRI Safety Information

Non-clinical testing demonstrated that the weighted CORFLO® Enteral Feeding Tubes are MR Conditional. A patient with this device can be scanned safely immediately after placement under the following conditions:

- Static magnetic field of 3-Tesla or less
- Spatial gradient magnetic field of 720-Gauss/cm or less
- MR system reported whole body averaged SAR of 3.0-W/kg (i.e., associated with a calorimetry measured value of 2.8- W/kg)

Under the conditions described above using a transmit/receive RF body coil the CORFLO® Enteral Feeding Tube produced the highest temperature change of +1.5 °C during MRI performed for 15-min in the 3-Tesla (3-Tesla/ 128-MHz, Excite, Software G3.0-052B, General Electric Healthcare, Milwaukee, WI) MR system.

Artifact Information: MR image quality may be compromised if the area of interest is relatively close to the position of the CORFLO® Enteral Feeding Tube. Therefore, optimization of MR imaging parameters to compensate for the presence of this device may be necessary.

Pulse Sequence	T1-SE	T1-SE	GRE	GRE
Signal Void Size	764-mm ²	89-mm ²	1,405-mm ²	684-mm ²
Plane Orientation	Parallel	Perpendicular	Parallel	Perpendicular

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COLOR YELLOW* is a Registered Trademark or Trademark of Avanos Medical, Inc.,
*Registered Trademark or Trademark of Avanos Medical, Inc., or its affiliates.
© 2018 AVNS. All rights reserved. 2024-12-03
Distributed in the USA by Avanos Medical Sales, LLC,
Avanos Medical, Inc., 5405 Windward Parkway, Alpharetta, GA 30004 USA.
Sponsored in Australia by Avanos Medical Australia Pty Ltd,