



Instructions of Use

Easy Access Sling

REPLACEMENT SLING RANGE

Reference: **EAS-IOU**

Product Codes: **M921532**


Revision: **21.11.2023.**

Copyright & Design Rights: **Silvalea.**



WARNING


Do **NOT** use clip sling attachments on a loop type spreader bar and Do **NOT** use loop sling attachments on a stud type spreader bar.

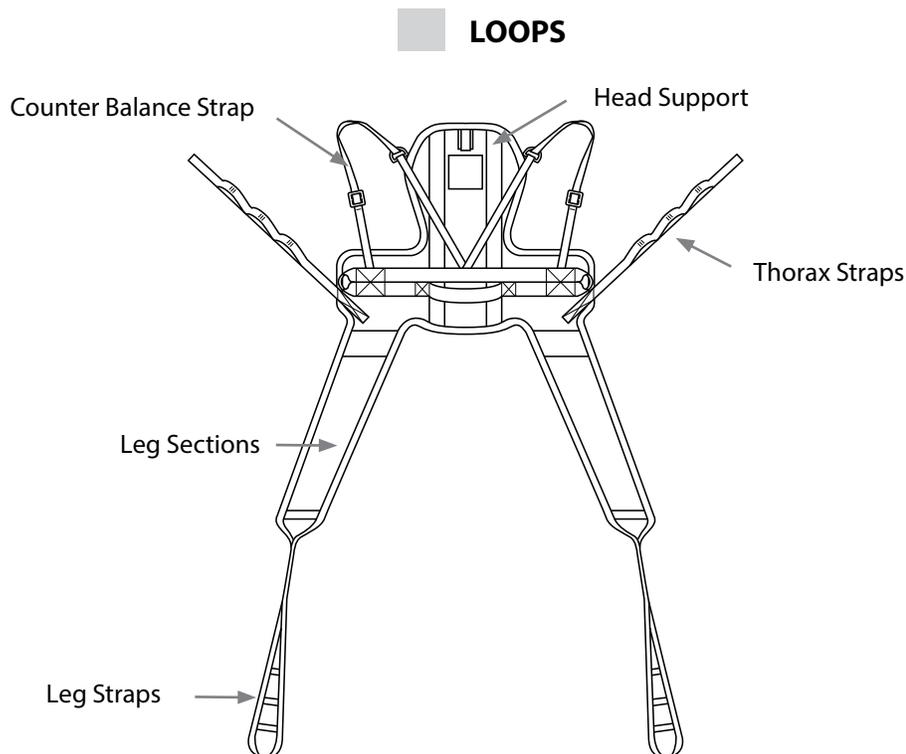

WARNING: Only use with suitable LOOP lift:




Sling Application:

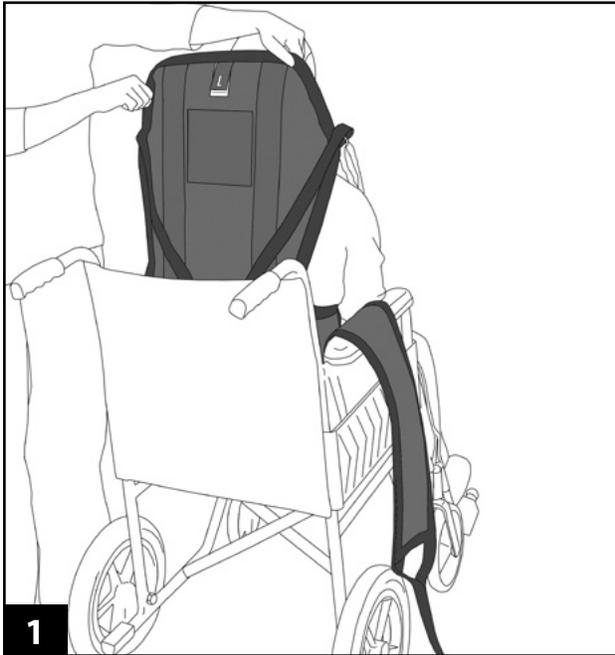


1 With the label facing outwards place the sling down behind the patient to approximately waist level.



2 Place the thoracic section around the patients trunk and cross over each thoracic strap and feed them through the loop on the opposite side.

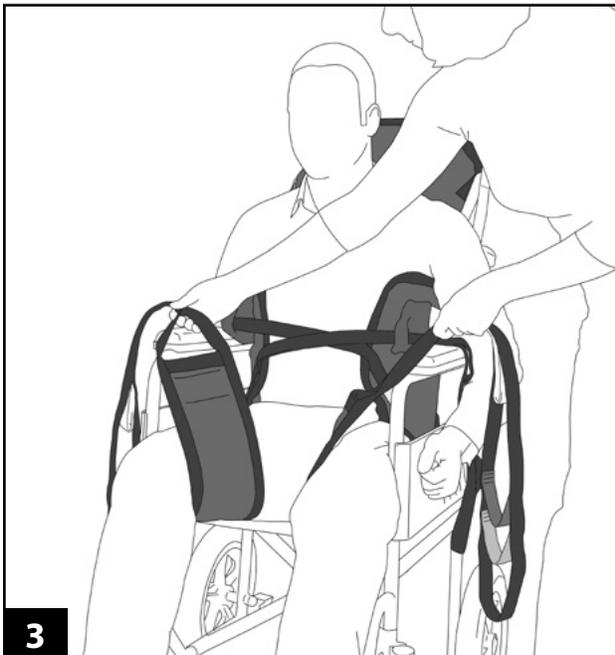


3 Glide each leg section under the patients thigh and bring between the legs. Feed one leg strap through the other to prevent the legs from abducting.

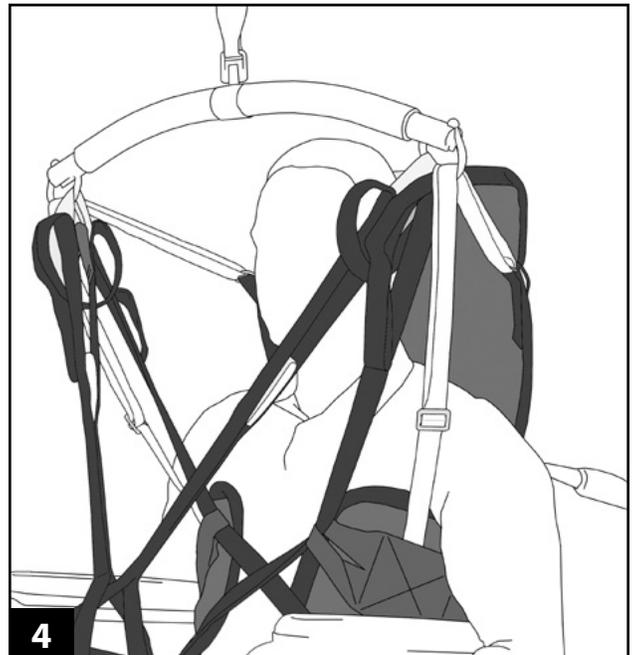


4 Attach both the legs and thoracic straps to the lifter using the appropriate color coded loops. With the arms on the outside (through the counterbalance strap) attach the strap to the lifter.

NB: Ensure the Counterbalance strap is applied last.

Sling Application Continued:

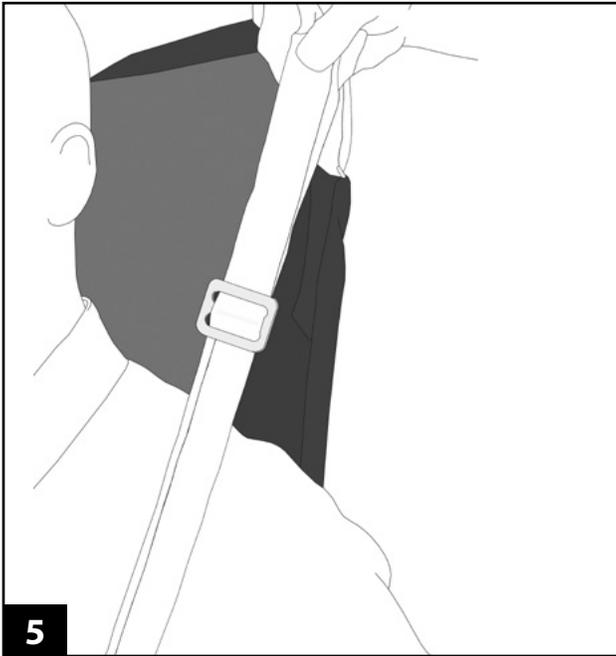

The counter balance strap can be easily adjusted to obtain the optimum position required by using the slider. There are 2 colored bartacks to help with positioning.

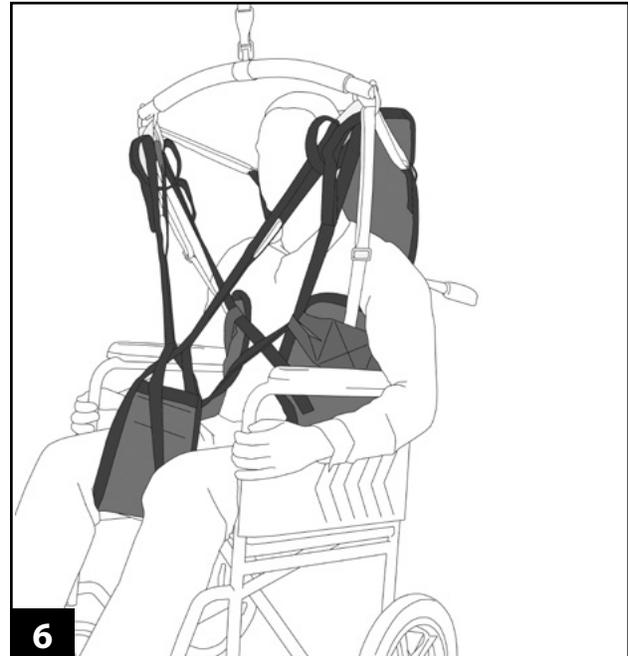

Once the sling has been fitted correctly to the patient and all straps are securely attached to the lifter raise the patient slightly and ensure that the thoracic strap is not too tight. Adjust if required.



You can now begin the transfer.



General Instructions: WARNING

Do not attempt to transfer the patient until all appropriate sling attachments are secured to the lift spreader bar and you are satisfied that the sling is the correct style, size and fit for the patient.

Correct sizing of a sling is important as a sling that is too large could cause the patient to slip through the aperture.

While every attempt is made to instruct a safe move you must conform to the current moving and handling procedures/regulations. If in any doubt, consult your Moving and Handling Advisor.

Service Label & Sling Checks:

All products shall be checked each time they are used for any signs of deterioration.

The service label is a visual guide to the product's check. Please conduct a thorough sling check in line with current safety requirements in your country. If in doubt consult your distributor for current sling checking protocols.

A risk assessment must be carried out to ensure that the correct size and style of sling has been met before use.

Pre-use sling check list:

- braiding/straps for deterioration and fraying
- material for holes/snagging
- compatibility with the lift
- safe working load of both the lift and sling/product (if applicable)
- label is legible

Remove from use immediately if there are any signs of deterioration to the sling or straps or the label is illegible.

Washing Guidelines:

- Follow laundry temperatures for washing and drying detailed on the product label. Please note: temperatures advised are MAXIMUM washing temperatures.
- Always ensure that all Hook & Loop fastenings are fully closed before washing.
- All products shall be laundered using non- biological detergents.
- DO NOT IRON

Product Care:

- Ensure that any Hook & Loop fastenings are closed when not in use. This will prevent unnecessary snagging of the sling material and prevent clogging of the Hook & Loop.
- Launder regularly.
- Store products in a suitable dry clean place.
- Silvalea will accept no liability for any additions/ adaptations made to our products that have NOT been carried out by ourselves. Doing so may cause serious injury or death to the patient.

Drying Guidelines:

- Follow the recommended drying temperatures as stated on product label. Please note: temperatures advised are MAXIMUM drying temperatures.
- Overheating a textile product can cause deterioration in the material components.
- Do not place the product on a convector heater or steam pipes.



SWL 272kg/600lb/42st



Silvalea Ltd a division of **Savaria**

