


ARJO

CONCERTO
SHOWER TROLLEY



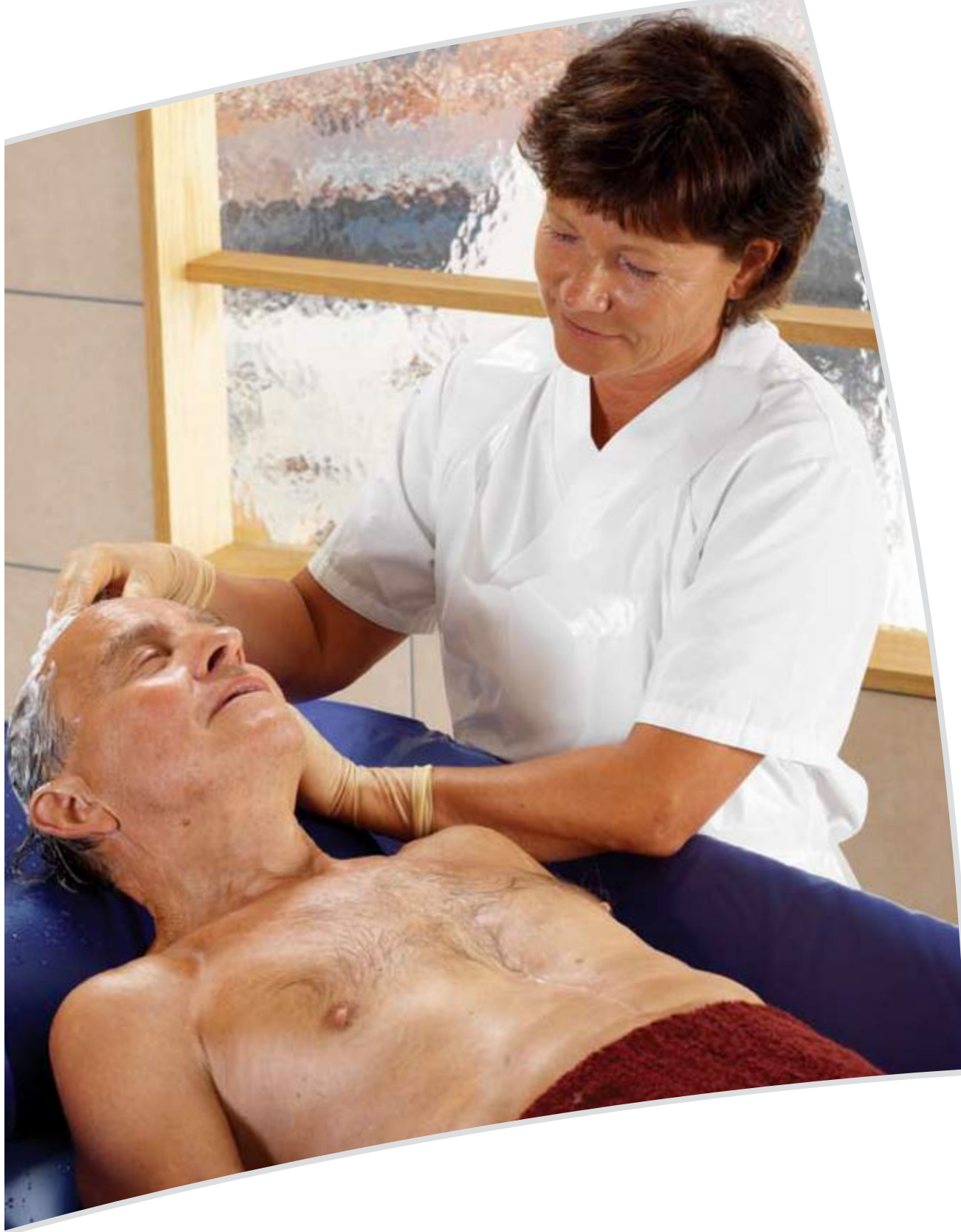

ASSISTED SHOWERING WITH EXTRA CARE

More and more healthcare facilities are turning to showering for residents' personal hygiene routines. *Concerto* offers a safe, well-proven system for showering dependent patients and residents in a range of healthcare environments.

Concerto can play an important role in supporting safe assisted showering and other hygiene routines at facilities for elderly care, acute care or special care.

This shower trolley is specifically designed as a safe platform to provide extra care for frail or disabled dependent patients or residents with no ability to support themselves, namely Emma, Edward, Emil and Elise in the Mobility Gallery.

Concerto is a proven showering solution that offers a far more hygienic and dignified alternative to a bed bath. Our shower trolley has been used successfully in healthcare facilities worldwide for many years to improve hygiene routines for dependent residents.

Patients are showered in a supine position on the comfortable and soft mattress. The trolley is stable and safe. Its fold down side supports have double safety latches to assure the security of residents and patients.

The wide, height-adjustable shower trolley offers ideal ergonomic working conditions, providing full body access for showering, hairwashing, pedicure and complete drying, as well as sufficient width to turn the individual easily.

Concerto is part of the comprehensive ARJO showering range. With *Concerto*, *Carendo* and *Carino*, ARJO provides healthcare facilities with optimum showering solutions for residents and patients at all levels of mobility.



Concerto is an ideal platform for ergonomic hygiene routines such as hairwashing.

ELDERLY CARE



The *Concerto* has a unique self leveling stretcher platform that allows overlapping on the bed to allow safe, ergonomic lateral transfers.



Smooth transport with a straight steering device for corridors.



The unique open chassis design on the *Concerto* allows easy safe transfer using sling lifts.



Regular foot care routines can be performed easily and efficiently.

Concerto improves the quality of life of dependent elderly residents by making showering routines safe, comfortable and dignified.

For many dependent residents in long-term care the only daily hygiene option is a bed bath, which is uncomfortable and undignified. A bed bath also puts a lot of strain on the caregiver's back and neck, as the procedure involves leaning over the bed.

Concerto provides an ergonomic solution for showering dependent residents. The height-adjustable function, operated via a handset on the battery-powered *Concerto* model, means caregivers can always use the shower trolley at the appropriate safe working height. For safe operating procedures, *Concerto* generally requires two caregivers.

Patient transfer to the shower trolley from the bed can be performed with a ceiling lift, mobile passive lifter or a horizontal transfer using *Maxi Slide* sliding sheets. Using the mobile *Concerto*, the resident can be transported to and from the shower area without any manual lifting.

The *Concerto* can be combined with flush sinks to allow easy, convenient splash-free showering by the caregivers. The stretcher slopes gently for good water run-off. Drain and run-off channels in the sealed mattress also speed water drainage. Used shower water is directed via a flexible drain hose to a floor drain or toilet, enabling a wider choice of rooms to be used as the shower area. After showering, the removable mattress and tiltable stretcher make cleaning/disinfection routines quick and easy.

Concerto is a proven solution for high-quality extra care in the assisted showering of dependent residents. Thousands of our shower trolleys are in daily use at elderly care facilities around the world.

In acute care, the *Concerto* can provide a good platform for cleaning wounds.

ACUTE CARE

At acute care facilities, *Concerto* can support routines that require stringent standards of hygiene.

Prior to surgical procedures, *Concerto* can act as an optimum station for the necessary pre-op washing for debridement and hygiene routines. For convenience, an optional IV bracket accessory is available. In wound or burn care, hygiene is a vital concern and *Concerto* offers a platform for good access when washing or treating affected parts of the patient's body.

A mobile solution can be of particular benefit in acute care. The *Concerto* can be moved to wherever it is needed in the facility. Our shower trolley, which can also serve as a useful nursing table, is available in three different lengths making it suitable for adults or children in acute care.



The soft mattress makes patients as comfortable as possible during pre-op hygiene procedures.



SPECIAL CARE

For dependent residents in special healthcare facilities, *Concerto* can offer a safe environment for assisted showering and a range of hygiene routines.

Special healthcare facilities cover a wide range of healthcare environments where residents may be dependent as a result of mental or physical disablement. With three different lengths available, *Concerto* enables assisted showering and hygiene routines for both children and adults.



Secure showering is essential for residents who have seizures. The extra-high side support option may be applicable for residents who are at risk.



Concerto allows staff to offer extra care by providing a safe platform for routines a resident cannot manage themselves, such as shaving.



ACCESSORIES

A range of *Concerto* accessories is available to maximize comfort and allow optimum adaptation for use in elderly, acute or special care.

The wedge pillow and adjustable back support both enhance patient comfort and give the caregiver better access for certain hygiene routines.

A resident or patient assessment may indicate that extra-high side supports will provide the right level of safety.

The IV bracket accessory is most applicable in acute healthcare environments.



Wedge pillow



Adjustable back support, three positions



Extra-high side supports



IV bracket accessory option

FEATURES



Electrical model, battery powered, comes complete with two batteries and a charger



Straight steering device



Flexible drain, drains to toilet or floor drain



Handset holder, electrical model only




Tiltable stretcher, allows for easy cleaning and storage




Drain hose holder

PRODUCT SPECIFICATIONS

Space requirement

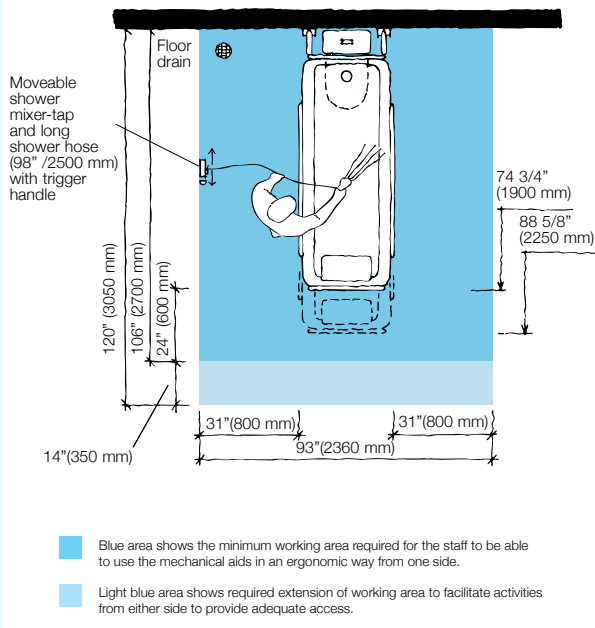

Mobility Gallery

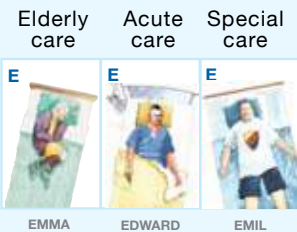

E, who is almost completely bedridden and totally dependant

Please contact ARJO for further information on the Mobility Gallery™.

Product Information

Lifting capacity	330 lbs (150 kg)
Weight	163 lbs (74 kg)
Weight (electric)	176 lbs (80 kg)
Width (stretcher)	25 5/8" (650 mm)
Total width (incl. side supports)	29 7/8" (760 mm)
Total length (different stretchers)	63" (1600 mm)
	74 7/8" (1900 mm)*
	88 1/2" (2250 mm)
Length of chassis (incl. wheel frame)	39 1/8" (995 mm)
Max internal length of chassis	32 1/4" (820 mm)
Max external width of chassis	27 1/2" (700 mm)
Min external width of chassis	11 7/8" (300 mm)
Height of chassis	5 7/8" (150 mm)
Height, top position floor to stretcher	35" (890 mm)
Height, top position floor to stretcher (electric)	34 1/4" (870 mm)
Height, low position floor to stretcher	22 1/2" (570 mm)
Stroke range	12 5/8" (320 mm)
Stroke range (electric)	11 3/4" (300 mm)
Color	Grey and blue
Equipped with head pillow and foldable side supports with safety lock	
Drain hose holders included	

Available as oil hydraulic (manual) version and battery powered (electrical) version.
*Adjustable +/- 2" (50 mm).

Use only ARJO Designed Parts, i.e. parts specifically designed for the purpose, on equipment and products supplied by ARJO. Our Conditions of Sale specify that ARJO is not liable for injuries attributable to the use of inadequate parts. Continuous development is part of our policy, and we therefore reserve the right to change specifications without notice.

Concerto, Carendo, Carino, Maxi Move, Maxi Slide and Mobility Gallery are trademarks belonging to ARJO International AG.

ARJO

...with people in mind.