

ACTA-EMBRACE™



Also available as a custom back. See page 52 for details.



FRONT VIEW



SHELL VIEW

FEATURES

Integrated trunk supports attached to a planar aluminum backrest shell may be adjusted for angle, height, and up to a total of 3" of width to create a custom supportive shape.

Aluminum stays on the backside of the two layer foam pad can be molded to support lumbar, scoliosis and kyphosis needs.

Aluminum stays are offset from center of the back, removing pressure points from the spine.

9" of total depth adjustment: 5" in front of canes and 4" behind the canes.

Comes with integrated chest harness strap guides.

Back comes with pre-drilled holes for optional extension strap guides.

Provides pelvic and scapular relief.

Compatible with head support pad: 2" wide x 1" tall hole pattern. Lightweight durable aluminum shell comes with pre-drilled holes for head support (Head Supports page 42).

Compatible with lateral trunk support hardware. A total of 11" of pad width adjustment (5.5" per side) and over 5" of height adjustment. Lateral trunk supports (E0956) are optional and available on page 41.

WEIGHT CAPACITY

10-18" Wide	250 lbs
20" Wide	300 lbs
22" Wide	350 lbs
24" Wide	400 lbs

STANDARD SIZES

	10" Wide	12" Wide	14" Wide	16" Wide	18" Wide	20" Wide	22" Wide	24" Wide
10" Height	10 x 10	10 x 12	10 x 14	10 x 16				
12" Height	12 x 10	12 x 12	12 x 14	12 x 16	12 x 18			
14" Height	14 x 10	14 x 12	14 x 14	14 x 16	14 x 18	14 x 20		
16" Height		16 x 12	16 x 14	16 x 16	16 x 18	16 x 20	16 x 22	
18" Height		18 x 12	18 x 14	18 x 16	18 x 18	18 x 20	18 x 22	18 x 24
20" Height			20 x 14	20 x 16	20 x 18	20 x 20	20 x 22	20 x 24
22" Height				22 x 16	22 x 18	22 x 20	22 x 22	22 x 24

*24" Wide backs require 4 Point Hardware (X & NX).

HARDWARE OPTIONS

C : Compass Quick Release 2 Point *See page 22 for details.


N : Compass Non-Quick Release 2 Point

X : Compass Quick Release 4 Point


NX : Compass Non-Quick Release 4 Point

COVER OPTIONS

STRETCH-AIR

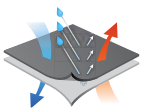

COMFORT-TEK

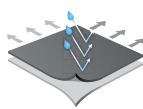

*See page 3 for details.



PART NUMBERS

Hardware	Stretch-Air	Comfort-Tek
COMPASS QUICK RELEASE 2-POINT	BK-EFXXSC-YY	BK-EFXXFC-YY
COMPASS NON-QUICK RELEASE 2-POINT	BK-EFXXSN-YY	BK-EFXXFN-YY
COMPASS QUICK RELEASE 4-POINT	BK-EFXXSX-YY	BK-EFXXFX-YY
COMPASS NON-QUICK RELEASE 4-POINT	BK-EFXXSNX-YY	BK-EFXXFNX-YY

XX & YY TO BE REPLACED. XX = HEIGHT & YY = WIDTH

STRAP GUIDE PART NUMBERS

Size	Comfort-Tek
EXTRA SMALL (UP TO 1.5" WIDE WEBBING)	CHG2-EB1
SMALL (UP TO 2.25" WIDE WEBBING)	CHG2-EB2