



ZM-310™

The Zippie ZM-310 combines our smallest, most reliable power wheelchair base with clinically proven pediatric seating and positioning options to provide children maximum accessibility and independence.



Cameron explores his world in the Zippie ZM-310

FEATURES AND BENEFITS

Our smallest mid-wheel-drive base is only 24" wide and 35" long with a turning radius of only 19", allowing your child to maneuver through narrow spaces

The ZM-310's reliable 4-pole, 6mph motors keep up with your child's most playful and energetic days




The ZM-310's Pediatric Specific Seating is available in sizes ranging from 10" X 10" to 16" X 19", and adjusts to accommodate growth

The ZM-310 is WC-19 Transit Approved for use as an occupant's seat in a motor vehicle when used with an approved Wheelchair Tiedown and Occupant Restraint System in accordance with SAE (J2249)

PG Drives Electronics and Optional Specialty Input Devices, including switch controls, head array, chin controls, and sip and puff, help make independent mobility possible



9" or 12" of Power Seat Elevate helps your child independently access tall shelves, cabinets, and light switches and provides the opportunity to interact with their peers at eye level

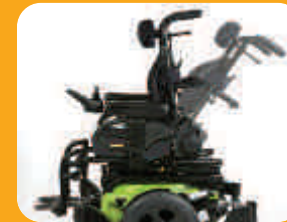

Manual Recline



Manual Tilt

50° of Power Tilt, 50° of Manual Tilt, and 40° of Manual Recline allow you or your child to change their position to help reduce the risk of skin breakdown, assist with feeding, digestion, respiratory function, and visual orientation, and decrease the potential of developing postural deviations

The ZM-310 is available with 15 shroud color options!



IMPORTANT NOTE: Print processes cannot accurately represent paint finishes. Actual colors may vary.

ZIPPIE ZM-310 SPECIFICATIONS

HCPCS Code	K0848, K0856	Seat Height	15.5" to 21.8"
Drive Wheel	Mid Wheel	Overall Base Width	24"
Max Speed	6 mph	Overall Length	35"
User Weight Capacity	165 lbs	Turning Radius	19"
Width	10" to 16"	Drive Wheels	13" with flat-free insert
Depth	10" to 19"	Transit Tested	WC-19



[Click here for printer-friendly version](#)