



HIPPOCAMPE®

MARATHON

Marathon and Trail



RUNNING ACCESSIBLE TO ALL

The Hippocampe Marathon and Trail allows you to practice running for leisure or in competition, and to **run in pairs on varied terrain.**

Lightweight spoked wheels make it easier to pick up speed, while the handlebar design and side grips optimise pushing the wheelchair. Extractable disc brakes provide effective braking.

Developed in close collaboration with runners, it has been designed to be as close as possible to your expectations.

The Hippocampe Marathon and Trail is the partner of your walks, hikes and races.



HIPPOCAMPE®

Marathon and Trail



TRAINING

Enjoy beautiful family outings on the road, walking or running. The slim spoked wheels make it easy to pick up speed.



MARATHON

10, 21 or 42km, challenge yourself as a duo or in team and share together the atmosphere of a run.



WALK

With the trail wheels, go for hikes in the forest or for walks in the countryside. The options allow you to adapt the Hiplocampe for a more or less sporty outing.



TRAIL RUN

The Hiplocampe Marathon and Trail allows you to participate in trails on seaside paths or in the forest thanks to the optional trail wheels.



ACTIVITIES

OPTIONS

HIPPOCAMPE®

MARATHON
Marathon and Trail



The Hippocampe Marathon and Trail is black and red and is available in 3 sizes: M, L and XL.

This wheelchair is offered with a range of options, which allow you to configure the chair to your needs.

CONFIGURATION ON ALL WHEELCHAIRS




Pair of armrests

The armrests, covered with foam, can be swivelled and retracted to facilitate transfers in the Hippocampe wheelchair.



Armrests nets

The armrest nets fit between the seat and armrests, to protect the user's hands from the wheel spokes.

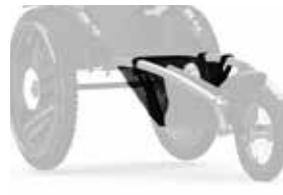

Footwell infill net

The footrest net fits between the seat and the footrest, to protect the user's feet.



Bottle and stand

The bottle is fixed thanks to its support on the push bar.

SUPPORT AND SAFETY OPTIONS



Mixed harness

Made of comfortable and water-resistant elastic neoprene, the 4-point harness facilitates the user's support and positioning.




Waist belt

Wide and comfortable, the waist belt keeps the user in a secure position at the bottom of the seat. A comfort foam protects the user from contact with the buckle.



Comfortable Plus Duo

The Stabilo neoprene positioning cushion can be shaped with a vacuum and provides additional support to the user by perfectly adapting to the body's shape and needs. It is waterproof.

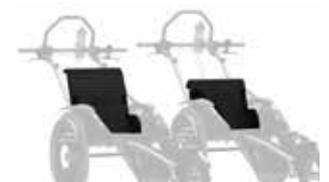

Height adjustable and reclining backrest

The height adjustable and reclining seat back frame makes it easy to vary the total height of the backrest according to the user's needs.



Adjustable headrest

The headrest is covered with comfortable foam and adapts to the user's position. It is adjustable forwards, backwards and in height.

COMFORT OPTIONS



Thermal cover

Comfortable and resistant, the thermal cover protects the user from cold and rain. Available in two sizes: M and L.



Rain canopy and legs protection

Suitable for rain, sun and wind, the rain canopy and legs protection is made of anti-UV materials. Thanks to its wide, transparent windows, the user can enjoy a panoramic view.



Umbrella

Compatible with the Hippocampe Beach and All-Terrain and the Hippocampe Marathon and Trail, the asymmetrical umbrella is wind resistant and protects from rain and sunrays.

LEISURE OPTIONS



Pair 24" trail wheels

The Hippocampe Marathon and Trail has been designed to be used on roads and, thanks to its optional trail wheels, on paths and trails.



Pushing horns

Installed on the handlebar, the bar-ends adaptation allows you to clip the handles of the chair vertically and position them according to the pusher's preferences.

BAGS



Transport bag

Useful when travelling, this bag allows you to transport the Hippocampe wheelchair during car or plane trips.



Backrest bag

The backrest bag, fitted on the back of the seat, allows to carry your personal belongings during your leisure activities and walks.



Bag for the 24 inches wheels

With this bag, trail or 24" road wheels are protected. Ideal for transport, this bag is equipped with two handles.



TECHNICAL INFORMATION

HIPPOCAMPE®

MARATHON® Marathon and Trail



UNIQUE COLOURS



Red and black

SIZE GUIDE



Size M: from 4' 7" - 1m40 to 5' 5" - 1m65

Size L: from 5' 5" - 1m65 to 6' 1" - 1m85

Size XL: from 6' 1" - 1m85

CE, FDA AND ISO 13485

VIPAMAT ensures a level of quality and performance that meets the certification requirements CE, of the FDA and ISO 13485.

TECHNICAL DATA


Seat: comfortable hydrophobic foam (3 cm thick) covered with a black fabric which dries quickly and prevents skin problems. Fabric and foam complying with the NFEN1021-2 fire resistance standard.

Frame: stainless steel and anodised aluminium covered with a foam cladding, comfortable to touch even when exposed to heat or cold.

Push bar: anodized aluminum tube with a handle covered with foam improving the comfort of the attendant.


Wheels: thin wheels for better performance on the road, and notched trail wheels for better grip on difficult terrains.

DIMENSIONS

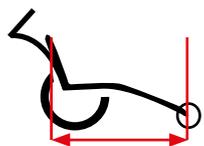

Seat length - feet

M: 32" - 82 cm
L: 36" - 92 cm
XL: 41" - 105 cm

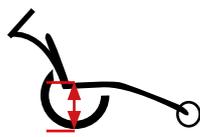


Seat height - ground

20 inches wheels: 13,4" - 34 cm
Marathon wheels 24 inches: 15.7" - 40 cm
Trail wheels 24 inches: 17.3" - 44 cm



Seat width

15.7" - 40 cm

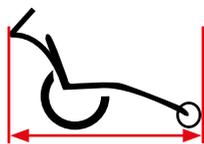

Total length

M: 69" - 175 cm
L: 73" - 185 cm
XL: 78" - 198 cm



Interior width (without armrests)

22" - 56 cm



Back seat height

14.5" - 37 cm to 18.5" - 47cm




Weight

Weight without option: 17 kg/37 lb
Max user weight: 100 kg/220 lb



Total width with trail wheels

27.6" - 70 cm

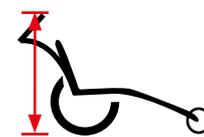

Total height

From 40.5" - 103 cm to 50.4" - 128 cm




Total length folded

M: 50.4" - 128 cm - L: 54.3" - 138 cm
XL: 59.4" - 151 cm



Total width with marathon wheels

27.6" - 70 cm

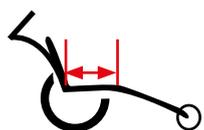

Seat depth

14.5" - 37 cm

